


ProInsert™

Simple

The ProInsert™ is a safe and easy-to-use device, designed for use with Nidacon products. The unique design, facilitates the safe and easy removal of the sperm pellet.

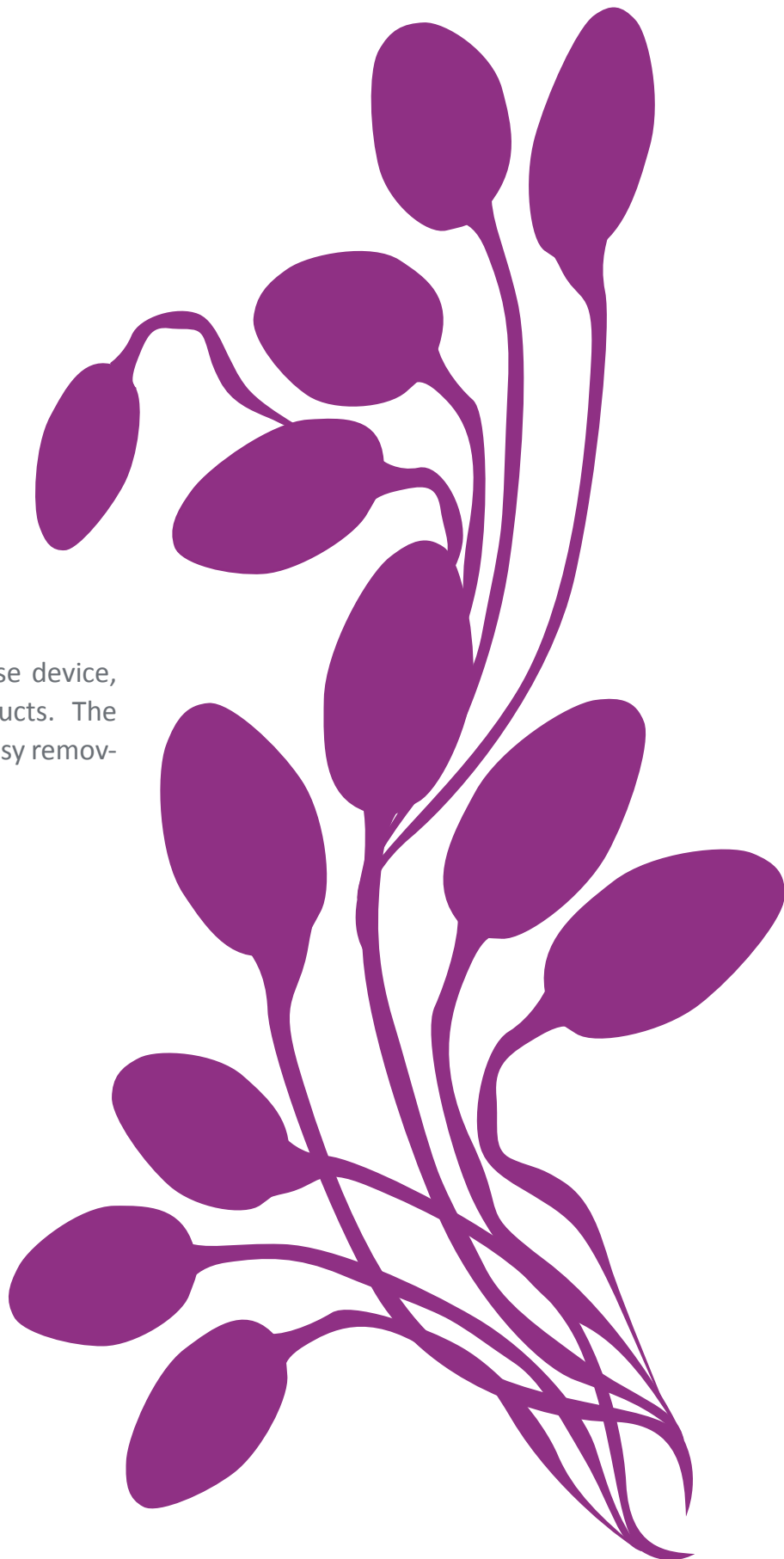

Makes the preparation of density gradient easy and fast

ProInsert™



Fast

ProInsert™ is used during preparation of human sperm from a semen sample using density gradient separation during assisted reproductive procedures. ProInsert™ facilitates density gradient preparation and pellet retrieval following density gradient separation.

Using the ProInsert™ reduced the time for completing a density gradient preparation of a semen sample, including sperm retrieval, by almost half.

Technical Information

ProInsert™ consists of an insert within a centrifuge tube, where the inside of the central channel never comes in contact with semen or the gradient.

Performance characteristics

Endotoxin level (EU/device)	<1.0
Human sperm survival 24 hours after density gradient preparation	>80% motility

The ProInsert™ has been sterilized by radiation.

Components

Centrifuge tube containing the ProInsert™
Pellet Retrieval Pipette

Storage conditions and shelf life

The unopened ProInsert™ Kit can be stored at room temperature until use.

24 months shelf life



Ordering Information

Volume

One pouch containing 5 individually packaged ProInsert™ kits

Article No.

PI15-5



www.nidacon.com

Nidacon International AB
E-mail contact@nidacon.com
Tel: +46 31 703 06 30 Fax: +46 31 40 54 15